

Yamato

Cold Trap

Model

CA300/800

Instruction Manual

- First Edition -

- Thank you for purchasing "Cold Trap, CA Series" of Yamato Scientific Co., Ltd.
- To use this unit properly, read this "Instruction Manual" thoroughly before using this unit. Keep this instruction manual around this unit for referring at anytime.



WARNING!:

Carefully read and thoroughly understand the important warning items described in this manual before using this unit.

Yamato Scientific Co. LTD.,


◆ Cautions in Using with Safety	1
• Explanation.....	1
• Table of Illustrated Symbols	2
• Fundamental Matters of “WARNING!” and “CAUTION!”	3
◆ Before Using this unit	5
• Requirements for Installation.....	5
◆ Description and Function of Each Part	8
• Main Unit and Control panel.....	8
◆ Installation Method	9
◆ Operation Method	10
• Procedure of Operation.....	10
• Cooling Curve (Reference Data)	11
◆ Handling Precautions	13
◆ Maintenance Method	14
• Daily Inspection and Maintenance	14
◆ Long storage and disposal	15
• When not using this unit for long term / When disposing	15
◆ In the Event of Failure	16
◆ After Service and Warranty	17
◆ Specification	18
◆ Wiring Diagram	19
◆ Replacement Parts Table	21
◆ Reference	22
• List of Dangerous Substances	22


Explanation

MEANING OF ILLUSTRATED SYMBOLS

Illustrated Symbols

Various symbols are used in this safety manual in order to use the unit without danger of injury and damage of the unit. A list of problems caused by ignoring the warnings and improper handling is divided as shown below. Be sure that you understand the warnings and cautions in this manual before operating the unit.

 **WARNING!** If the warning is ignored, there is the danger of a problem that may cause a serious accident or even fatality.

 **CAUTION!** If the caution is ignored, there is the danger of a problem that may cause injury/damage to property or the unit itself.

Meaning of Symbols



This symbol indicates items that urge the warning (including the caution). A detailed warning message is shown adjacent to the symbol.



This symbol indicates items that are strictly prohibited. A detailed message is shown adjacent to the symbol with specific actions not to perform.



This symbol indicates items that should be always performed. A detailed message with instructions is shown adjacent to the symbol.


Cautions in Using with Safety

Table of Illustrated Symbols

Warning



Warning,
generally

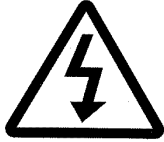

Warning,
high voltage



Warning,
high temperature



Warning,
drive train



Warning,
explosive

Caution



Caution,
generally



Caution,
electrical shock



Caution,
scald



Caution,
no road heating



Caution,
not to drench



Caution,
water only



Caution,
deadly poison

Prohibit



Prohibit,
generally



Prohibit,
inflammable



Prohibit,
to disassemble



Prohibit,
to touch

Compulsion



Compulsion,
generally



Compulsion,
connect to the
grounding
terminal



Compulsion,
install on a flat
surface



Compulsion,
disconnect the
power plug



Compulsion,
periodical
inspection

Cautions in Using with Safety

Fundamental Matters of “WARNING!” and “CAUTION!”

WARNING!

Do not use this unit in an area where there is flammable or explosive gas

Never use this unit in an area where there is flammable or explosive gas. This unit is not explosion-proof. An arc may be generated when the power switch is turned on or off, and fire/explosion may result. (Refer to page22 “List of Dangerous Substances”.)

Always ground this unit

Always ground this unit on the power equipment side in order to avoid electrical shock due to a power surge.

If a problem occurs

If smoke or strange odor should come out of this unit for some reason, turn off the power key right away, and then turn off the circuit breaker and the main power. Immediately contact a service technician for inspection. If this procedure is not followed, fire or electrical shock may result. Never perform repair work yourself, since it is dangerous and not recommended.

Do not use the power cord if it is bundled or tangled

Do not use the power cord if it is bundled or tangled. If it is used in this manner, it can overheat and fire may be caused.

Do not process, bend, wring, or stretch the power cord forcibly

Do not process, bend, wring, or stretch the power cord forcibly. Fire or electrical shock may result.

Pay special attention to the measure for flammability and handling of flammable solvent

Leaving at the temperature higher than the room temperature may vaporize the flammable material (ethanol, etc.). There might be the case that some flammable liquid might be vaporized at the temperature lower than the room temperature. The result of such careless handling could cause the fire or explosion. Do provide the vaporization with enough during the operation.

Do not disassemble or modify this unit

Do not disassemble or modify this unit. Fire or electrical shock or failure may be caused.

Cautions in Using with Safety

Fundamental Matters of “WARNING!” and “CAUTION!”

CAUTION!

During a thunder storm

During a thunderstorm, turn off the power key immediately, then turn off the circuit breaker and the main power. If this procedure is not followed, fire or electrical shock may be caused.

Do not touch the condensed liquid in the cooling coil and trap bath

Since the condensed liquid in the cooling coil and trap bath is stayed with low temperature, never touch it so as to preventing from getting frostbite on your hands.

Do not touch the cooling fin with bare hands

Do not touch the cooling fin with bare hands during maintenance, for the edge of the cooling fin is too sharp to cut your hand.

Requirements for Installation

WARNING!

1. Always ground this unit



- Connect the power plug to a receptacle with grounding connectors.



- Do not forget to ground this unit, to protect you and the unit from electrical shock in case of power surge. Choose a receptacle with grounding connectors as often as possible.
- Do not connect the grounding wire to a gas pipe, or by means of a lightning rod or telephone line. A fire or electrical shock will occur.


2. Choose a proper place for installation



- Do not install this unit in a place where:
 - ◆ Rough or dirty surface.
 - ◆ Flammable gas or corrosive gas is generated.
 - ◆ Ambient temperature below 5°C or above 35°C. (CA300)
 - ◆ Ambient temperature below 5°C or above 30°C. (CA800)
 - ◆ Ambient temperature fluctuates violently.
 - ◆ There is direct sunlight.
 - ◆ There is excessive humidity and dust.
 - ◆ There is a constant vibration.




- Install this unit on a stable place with the space as shown below.

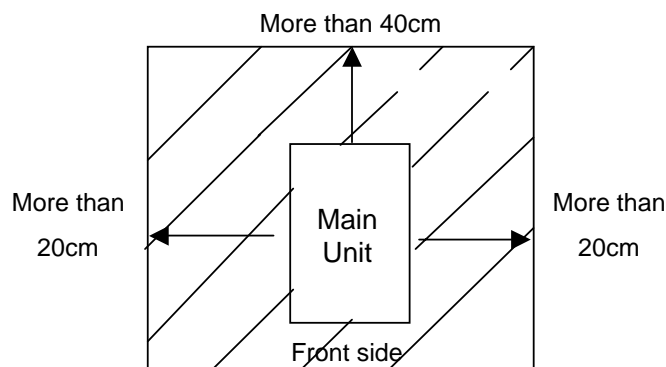

3. Do not use this unit in an area where there is flammable or explosive gas



- Never use this unit in an area where there is flammable or explosive gas. This unit is not explosion-proof. An arc may be generated when the power switch is turned ON or OFF, and fire/explosion may result.
- To know about flammable or explosive gas refer to page22 "List of Dangerous Substances".

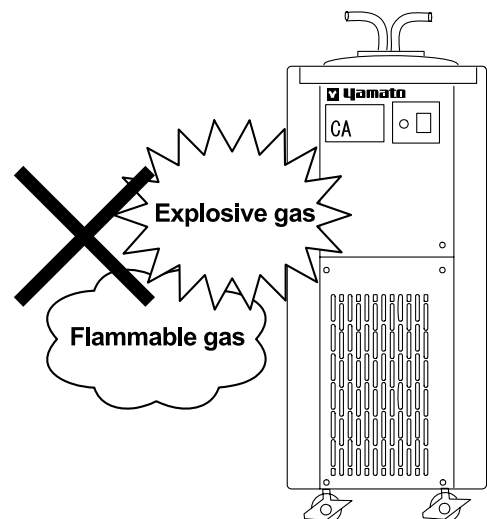


Requirements for Installation

4. Do not modify




- Modification of this unit is strictly prohibited. This could cause a failure.

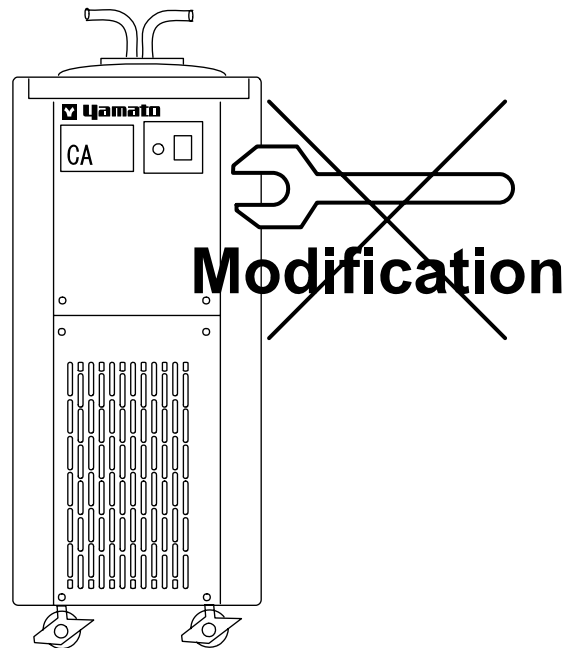

5. Do not topple or tilt this unit



- Set this unit to the flattest place. This unit incorporates the refrigerator. Do not topple or tilt it.

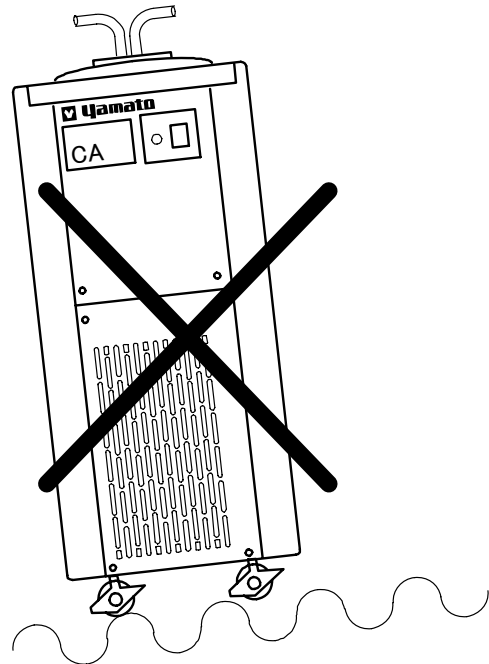

Requirements for Installation

CAUTION!

6. Use specified receptacle for power source



- Choose a correct power distribution board or receptacle that meets the unit's rated electric capacity.

Electric capacity: CA300: 100V AC, 50/60Hz, 6A
CA800: 100V AC, 50/60Hz, 10A

NOTE)

Starburst connection with a branching receptacle or extended wiring with a cord reel lowers electrical power voltage, which may cause the degradation of refrigeration capability.

7. Before/after installing



- It may cause injure to a person if this unit falls down or moves by the earthquake and the impact. etc..To prevent, take measures that the unit cannot fall down, and not install to busy place.
- Though this unit has the air-cooled refrigerator, the device exhausts the heat. Do provide the vaporization with enough so as not to raise the ambient temperature caused by the exhaust of the heat, or install this unit with its air controlled completely. If the ambient temperature becomes high, the operation efficiency becomes worse, and could cause the malfunction of the device by high temperature and humidity.

8. Handling of power code



- Do not entangle the power cord. This will cause overheating and possibly a fire.
- Do not bend or twist the power cord, or apply excessive tension to it. This may cause a fire and electrical shock.
- Do not lay the power cord under a desk or chair, and do not allow it to be pinched in order to prevent it from being damaged and to avoid a fire or electrical shock.
- Keep the power cord away from any heating equipment such as a room heater. The cord's insulation may melt and cause a fire or electrical shock.




- If the power cord becomes damaged (wiring exposed, breakage, etc.), immediately turn off the power at the rear of this unit and shut off the main supply power. Then contact your nearest dealer for replacement of the power cord. Leaving it may cause a fire or electrical shock.
- Connect the power plug to the receptacle which is supplied appropriate power and voltage.


Description and Function of Each Part

Main Unit and Control panel

Main unit

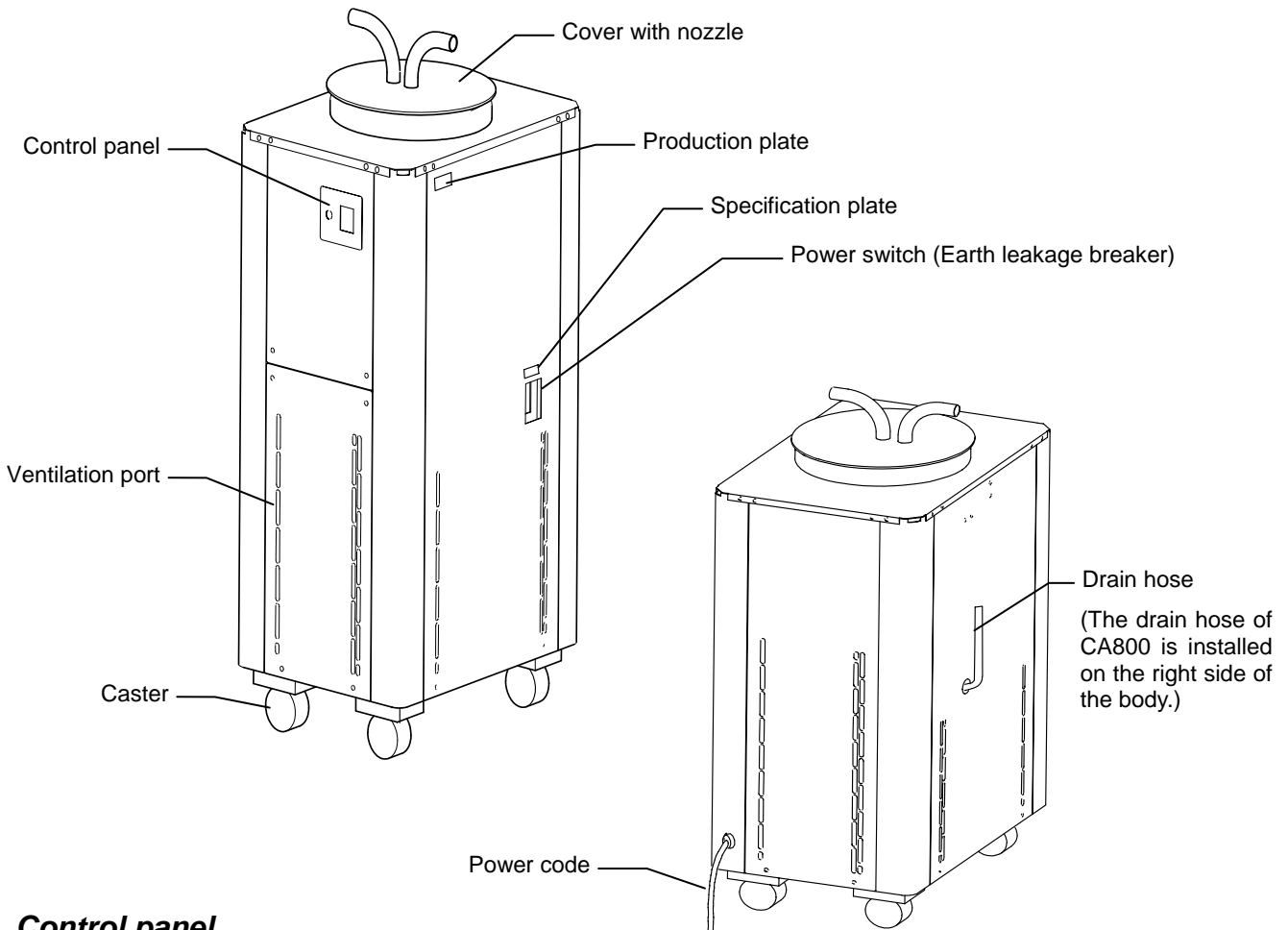



Control panel

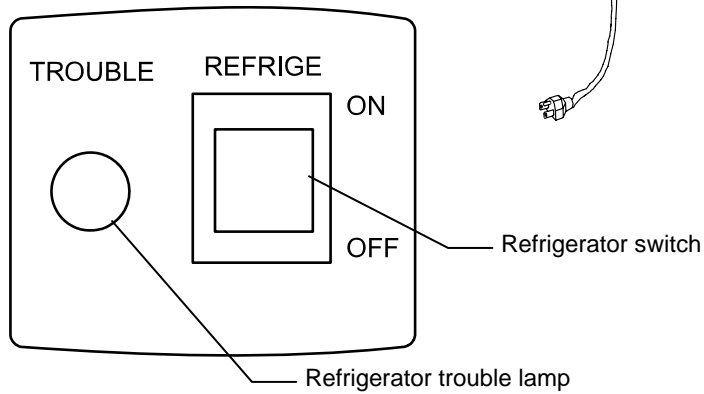

Function of each part

Part Name	Function
Earth Leakage Breaker :	This is the power switch Turns ON/OFF the main power.
Refrigerator Switch :	The refrigerator activates by pressing "ON" of this switch.
Refrigerator Trouble Lamp :	This lamp lights on when the refrigerator is in overload.

1. Unlock the stopper of the caster.
Pulling up the lever of the stopper for caster releases the lock.
(Only the two casters in front of the unit are attached the stopped.)
2. Move the device to the place to be installed.
 - ❖ If there is a step on the floor, the too strong impact is given to the caster, and could give the damage. In that case, move the device by lifting at the step.
3. When the installation place is determined, pull down the lever of the stopper for caster, and lock them.
4. Drain Tap Check
 - Check whether the drain tap is detached or not.




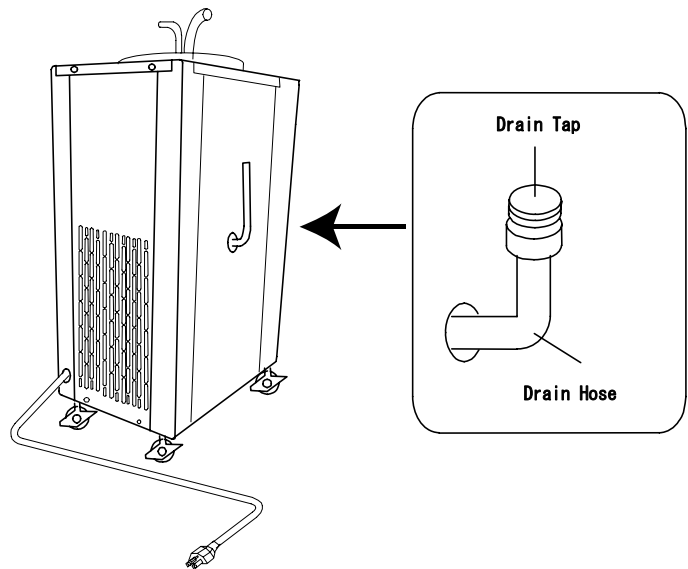
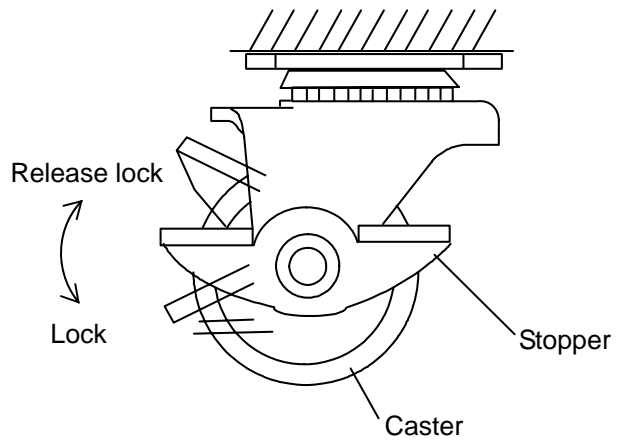


5. Power Plug Connection
 - Check the power of the earth leakage breaker is turned "OFF", and plug the power cord in the receptacle.

Procedure of Operation

Read through the following procedure if the liquid to be trapped is water type.


Use the glass condenser set (optional accessory) for the acid or organic solvents. Refer to the operation manual attached to the glass condenser for using this glass condenser.
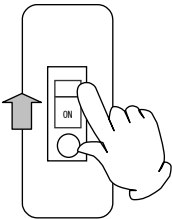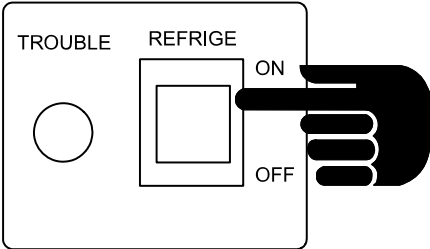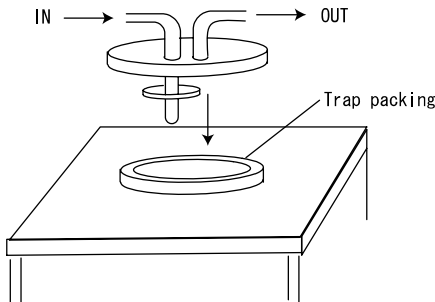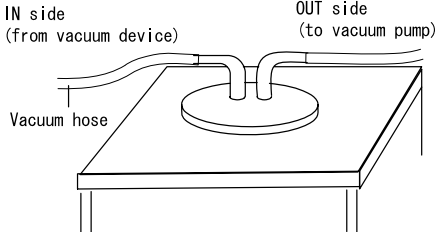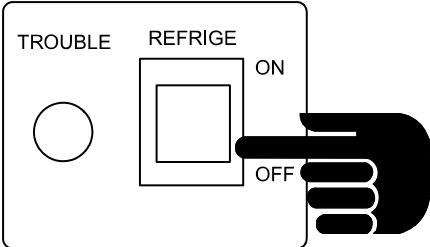
<p>1 Supply the power</p> 	<p>1. Turn "ON" the earth leakage breaker.</p>
<p>Operation of the refrigerator</p> 	<p>2. Turn "ON" the refrigerator switch. The refrigerator activates.</p> <ul style="list-style-type: none"> ❖ There is the case that the starting sound of the refrigerator might become higher depending on the operation initial status and ambient temperature. However, this phenomenon is not abnormal.
<p>2 Setting of cover with nozzle</p> 	<p>1. Check whether there is a dent or dust on the trap packing or cover with nozzle or not.</p> <ul style="list-style-type: none"> ❖ If any dent or dust is detected, the vacuum degree of the inside of the trap bath becomes worse. <p>2. Mount the cover with nozzle on the trap packing as checking the "IN" and "OUT" of the nozzle direction.</p> <ul style="list-style-type: none"> ❖ After setting the cover with nozzle, check the vacuum degree of the sealing parts. Insufficient sealing could cause the air leakage from the trap packing.
<p>3 Vacuum hose connection</p> 	<p>Connect the cover with nozzle with its "IN" side connected to vacuum equipment and with its "OUT" side connected to the vacuum pump using the vacuum hose.</p> <ul style="list-style-type: none"> ❖ There is "IN" and "OUT" side for the cover with nozzle. Do not connect to inverse side. ❖ The vacuum hose is not an attached accessory. Prepare it by yourself. Apply the vacuum hose with the inner diameter 12mm or 15mm.
<p>4 Operation of the vacuum pump</p>	<p>Turn "ON" the vacuum pump. The vacuum pump collects the steam from the vacuum equipment to the trap bath.</p>
<p>5 End of operation</p> 	<p>1. Stop operating the device.</p> <p>2. Remove the vacuum hose from the cover with nozzle.</p> <p>3. Turn "OFF" the vacuum pump.</p> <p>4. Melt the ice in the trap leaving in the natural temperature or by pouring the hot water. After melted all ice, drain all the water from the device by removing the drain tap.</p> <p>NOTE) Be sure to turn "OFF" the vacuum pump switch after turning back the pressure of the vacuum pump to normal one for preventing from the inverse flow of the oil.</p>


Cooling Curve (Reference Data)

CA300

■ Ethanol 3liters

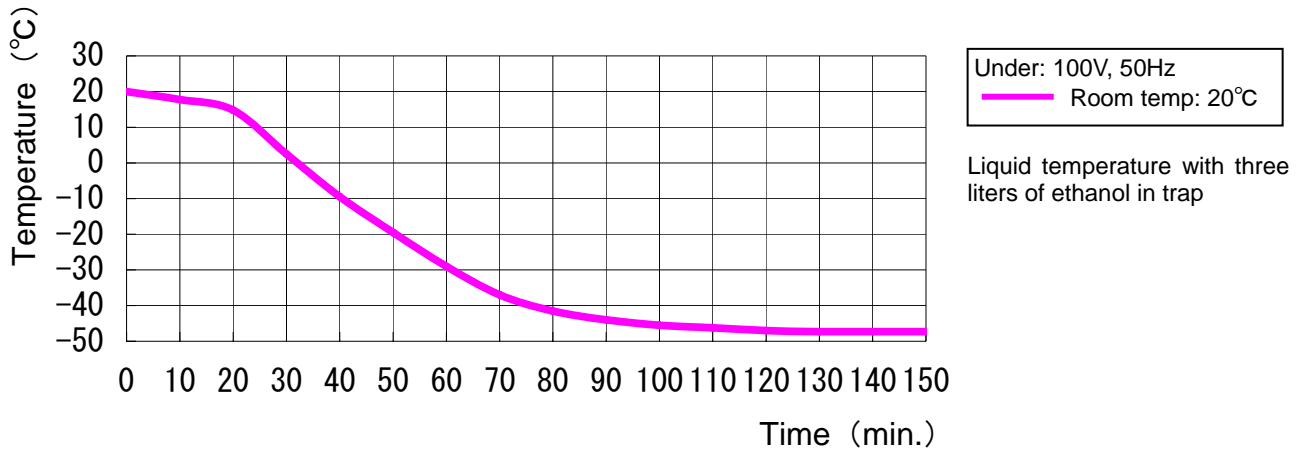


■ Cooling part of the trap (surface)

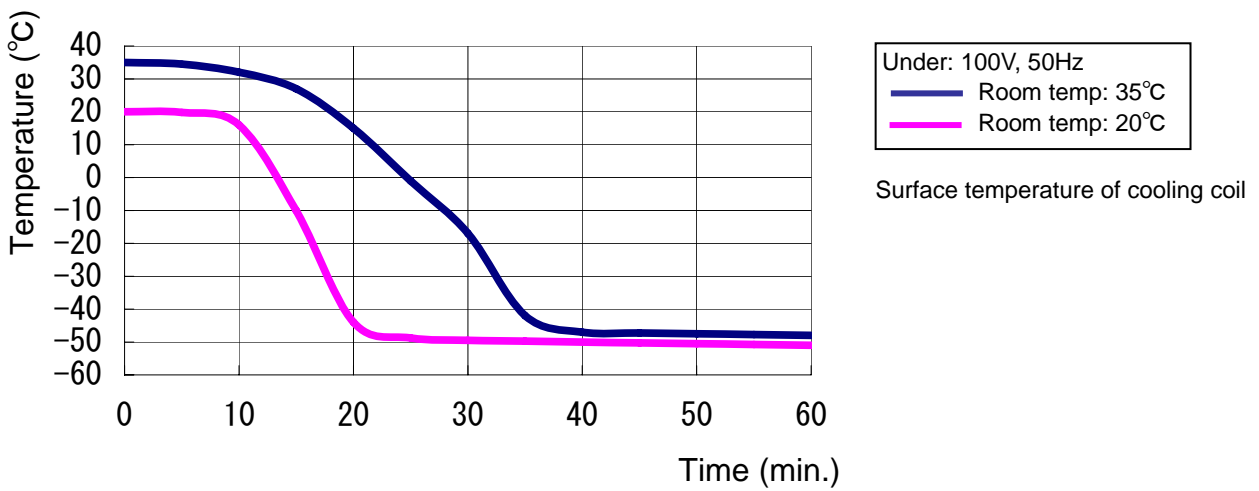


Cooling Curve (Reference Data)

CA800

■ Ethanol 3liters

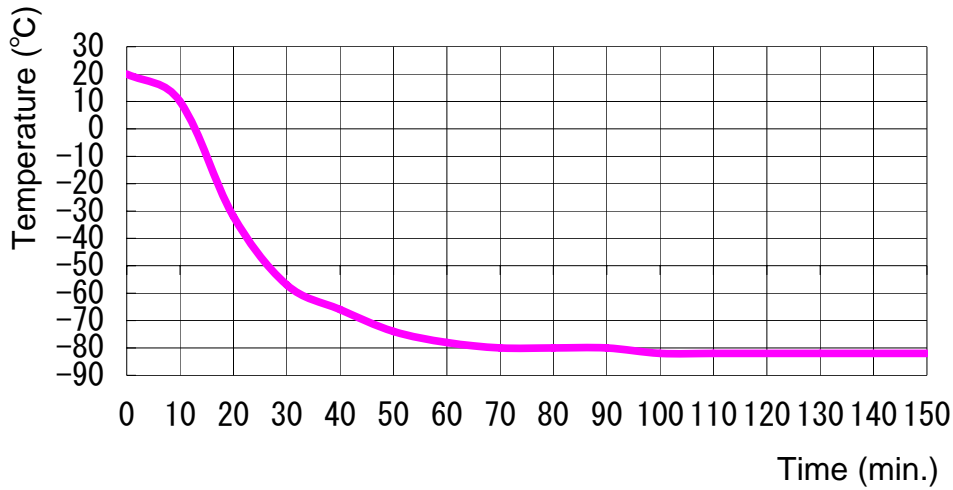


■ Cooling part of the trap (surface)

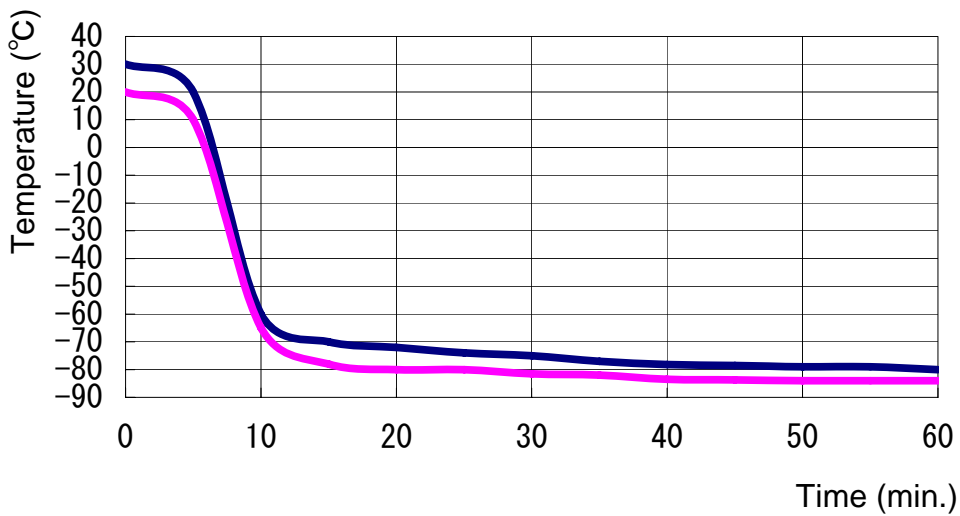


WARNING!

If a problem occurs


-  If smoke or strange odor should come out of this unit for some reason, turn off the power key right away, and then turn off the breaker and the main power. Immediately contact a service technician for inspection. If this procedure is not followed, fire or electrical shock may result. Never perform repair work yourself, since it is dangerous and not recommended.

Measure for flammability and handling of flammable solvent


-  This unit is not designed as the explosion-proof construction. Pay special attention to the handling of the sample to be handled with this unit on the consumption with the explosive material, flammable material, and similar ones. The flammable material may be vaporized by leaving it at the temperature higher than room temperature, and could cause the fire or explosion. When handling such material, provide ventilation with enough before the operation. (Refer to page22 "List of Dangerous Substances".)

CAUTION!


Maximum trap capacity

-  The max. trap capacity for the water type liquid is around 0.9kg (CA300) or 0.8kg (CA800). Inserting too much trap to this unit could cause the clogging of "IN", "OUT" nozzle.


Do not step on this unit

-  Do not step on this unit. It will cause injury if this unit fall down or break.


Do not put anything on this unit

-  Do not put anything on this unit. It will cause injury if fall.

During a thunder storm

-  During a thunderstorm, turn off the power key immediately, then turn off the circuit breaker and the main power. If this procedure is not followed, fire or electrical shock may be caused.

Countermeasure for stop operation during night or long-term stop

-  In case of stopping operation during night or long-term, toggle the power switch to "OFF".


Daily Inspection and Maintenance

WARNING!

- Disconnect the power cable from the power source when doing an inspection or maintenance unless needed.
- Perform the daily inspection and maintenance after returning the temperature of this unit to the normal one.
- Do not disassemble this unit.
- Do not touch the cooling fin with bare hands.


CAUTION!

- Use a well-drained soft cloth to wipe dirt on this unit. Do not use benzene, thinner or cleanser for wiping. Do not scrub this unit. Deformation, deterioration or color change may result in.

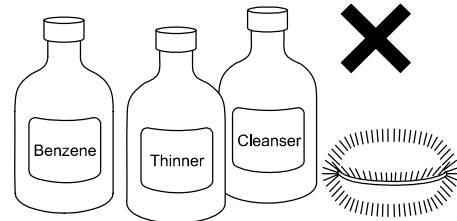

Monthly maintenance


- Check the earth leakage breaker function.
 1. Connect the power cord.
 2. Turn the breaker on.
 3. Push the red test switch by a ballpoint pen etc.
 4. If there is no problem, the earth leakage breaker will be turned off.

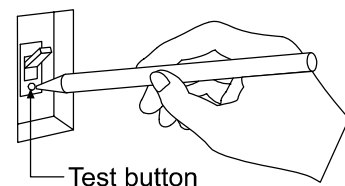

Cleaning of cooling fin

- Clogging of the cooling fin could cause the deterioration of the cooling performance, and also cause the malfunction of the refrigerator. The clogged status differs depending on the surrounding condition or operation time. Clean the cooling fin periodically.
- ❖ Loosen the mounting screws (4 screws) of the ventilation port cover, remove the cover of the ventilation port, and remove the dust attached to the surface of the cover using the vacuum cleaner.
- ❖ After cleaning the cooling fin, attach it in inverse procedure.

 Take care not to crush the fin during cleaning.

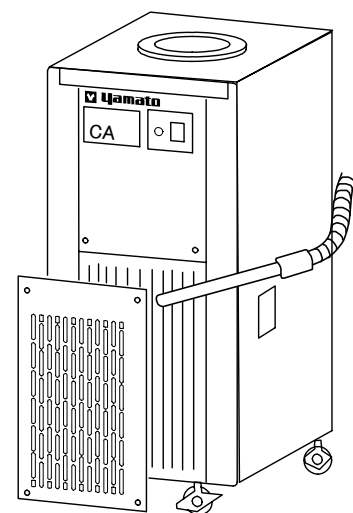

For any questions, contact the dealer who you purchased this unit from, or the nearest sales division in our company.

When not using this unit for long term / When disposing

CAUTION!

When not using this unit for long term...

- Turn off the power and disconnect the power cord.

WARNING!

When disposing...

- Keep out of reach of children.
- Treat as large trash.

Environmental protection should be considered

We request you to disassemble this unit as possible and recycle the reusable parts considering to the environmental protection. The feature components of this unit and materials used are listed below.

Component Name	Material	
Parts of Main Unit		
Casing	Bonderizing steel plate baked with melamine resin coating	
Inner bath, Cover	Stainless steel SUS304	
Production plates	Polyester (PET) resin film	
Corner	Alkylbenzenesulfied (ABS) resin	
Trap packing	Silicon rubber	
Electrical Parts		
Switches, Relays	Composites with resin and others	
Power cord & wiring materials and others	Composites with synthetic rubber, copper, nickel and others	
Parts of Refrigeration System		
Refrigerator unit	Iron, Synthetic oil and others	
Cooling fin	Aluminum, Copper	
Parts of Piping		
Drain hose	Natural rubber	
Drain tap	Polyacetal	
Drain tap holder	Bonderizing steel plate baked with melamine resin coating	
Pipe cover	Polyurethane sponge	
Pipe	Copper	
Sealed Cooling Medium for Refrigerator		
Cooling medium	R404A	Ask the specialist for the dealing of cooling medium.
	HFC23	

Trouble Shooting

Condition	Possible Causes
The device does not start when turning on the power switch.	<ul style="list-style-type: none"> • Power plug is not connected to the receptacle correctly. • Power failure. • Earth leakage breaker is turned to "OFF"
Not fallen the temperature.	<ul style="list-style-type: none"> • The cooling fin is clogged. • The cooling liquid is overheated. • The ambient temperature is exceeding 30°C(CA800) or 35°C (CA300). • The peripheral of the ventilating port is shut down.

In lighting on the lamp,

Error Sign	Cause	Remedy
Turned on the refrigerator trouble lamp	Overload of refrigerator	<ul style="list-style-type: none"> • Immediately turn off the power, remove the cause of the error referring to the "Trouble Shooting" (Not fallen the temperature) above, and turn on the power again after passing the certain time. • In case of lighting the alarm lamp again, make a call to the service office.

In the case if the error other than listed above occurred, turn off the power switch and primary power source immediately. Contact the shop of your purchase or nearest Yamato Scientific Service Office.

In Case of Request for Repair

If the failure occurs, stop the operation, turn OFF the power switch, and unplug the power plug. Please contact the sales agency that this unit was purchased, or the Yamato Scientific's sales office.

< Check following items before contact >

- ◆ Model Name of Product
 - ◆ Production Number
 - ◆ Purchase Date
 - ◆ About Trouble (in detail as possible)
- } See the production plate attached to this unit.

Minimum Retention Period of Performance Parts for Repair

The minimum retention period of performance parts for repair of this unit is 7 years after discontinuance of this unit.


The "performance part for repair" is the part that is required to maintain this unit.

Model		CA300	CA800
Method		Direct Trap or glass trap (Optional accessory)	
Performance	Dehumidifying capacity	Max. 0.9Kg (Water type liquid)	Max. 0.8Kg (Water type liquid)
	Lowest reached temperature	-45°C	-80°C
	Time required to reach lowest temperature	20min.	30min.
Configuration	Refrigerator	Air cooling, 300W	Air cooling, 350W × 2
	Cooling medium	R404a	R404a/R23
	Cooling coil	Inner diameter: 90mm, SUS304	
	Cover	SUS304, Outer diameter: 17.6mm, with nozzle	
	Bath material	SUS304, Cylindrical	
	Dimensions of bath	φ 153 × H235mm	
	Capacity	Approx. 4L (liquid quantity: 3L)	
Usable ambient temperature		5 to 35°C	5 to 30°C
Standards	External dimensions	W284 × D404 × H700 mm	W405 × D455 × H850 mm
	Weight	Approx. 32Kg	Approx. 54Kg
	Power supply	100V AC, 50/60Hz, 6A	100V AC, 50/60Hz, 10A
Accessory		Instruction manual	

Optional Accessories (common for CA300/800)


Name		Product code	Notes
Glass condenser set OCA10		2221487	Collection of acid and organic solvents
Tube with different dia. for rubber tube	Brass φ 30 × φ 18	242185	For joint of parts with different diameter
	Brass φ 18 × φ 12	242186	
	SUS φ 30 × φ 18	241496	
	SUS φ 18 × φ 12	241497	
Cover made of SUS		200000	Cover of trap bath

CA300

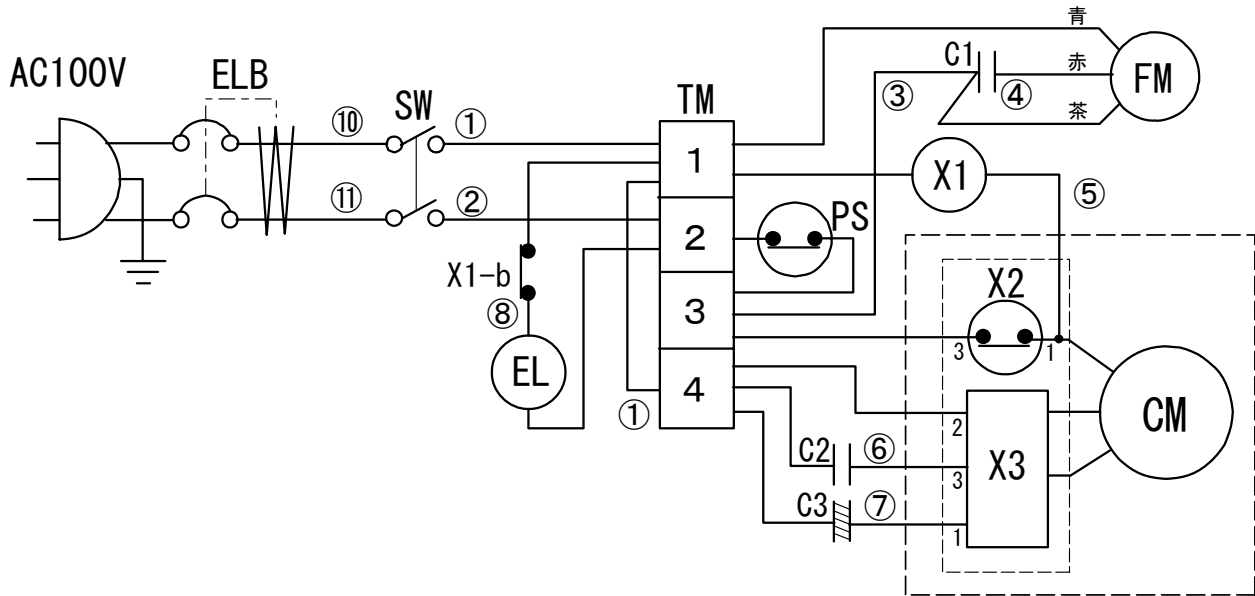

Symbol	Part name	Symbol	Part name
ELB	Earth leakage breaker	PS	Pressure switch
SW	Refrigerator switch	C1	Operation condenser (FM)
TM	Terminal block	C2	Operation condenser (CM)
EL	Trouble lamp	C3	Start condenser
X1	Relay	FM	Fan motor
X2	Overload relay	CM	Compressor motor
X3	PTC start relay		

CA800

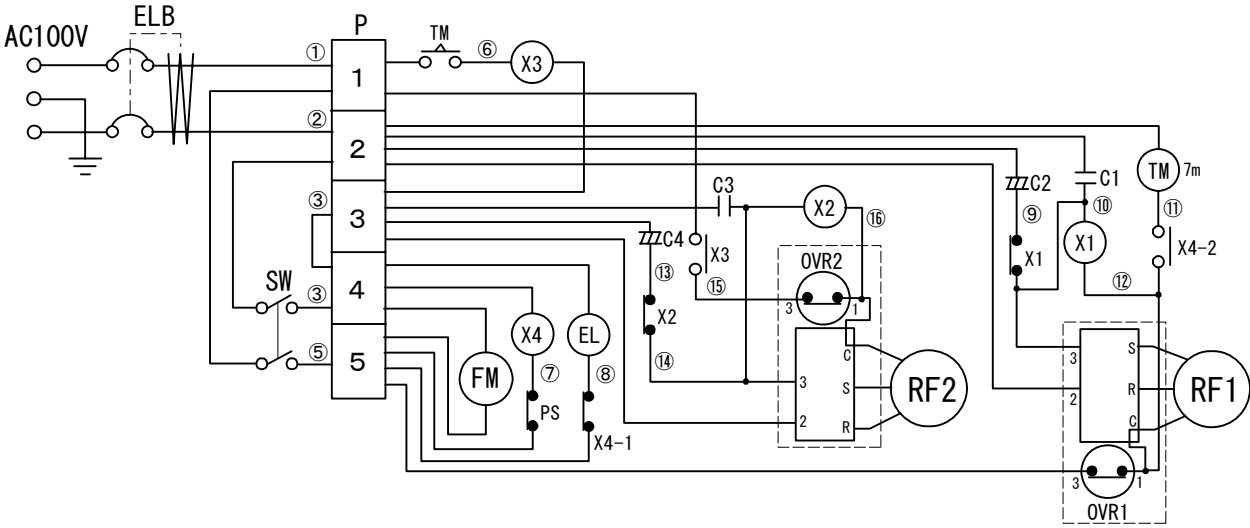

Symbol	Part name	Symbol	Part name
ELB	Earth leakage breaker	OVR1, OVR2	Overload relay
SW	Refrigerator switch	C1, C3	Operation condenser
TM	Terminal block	C2, C4	Start condenser
EL	Trouble lamp	TM	Time relay
X1, X2	Start relay	FM	Fan motor
X3, X4	Relay	RF1, RF2	Refrigerator
PS	Pressure switch		

Replacement Parts Table

CA300

Part Name	Code No.	Specification	Manufacturer
Packing	CA300-30150	Silicon	Yamato Scientific
Refrigerator unit	3-01-003-6005	100V/300W UF-NS300L-YL	Yamato Scientific
Capillary	3-07-002-0004	Inner dia. ϕ 0.8 \times outer dia. ϕ 1.8	Yamato Scientific
Refrigerator switch	2-01-027-0005	HLN215A	Fujisoku
Trouble lamp	2-09-006-0001	BN-5701	Satoh Parts
Earth leakage breaker	2-06-000-0019	FG32R/15-30MA15A	Fuji Denki
Relay	2-05-008-0003	AP3124K	Matsushita
Drain hose	3-04-001-6022	ϕ 6 \times ϕ 18	Yamato Scientific
Terminal block	LT00035672	MKH-250ABC-4P	Terminal
Drain tap	CA300-40200	Polyacetal	Yamato Scientific

CA800

Part Name	Code No.	Specification	Manufacturer
Packing	CA300-30150	Silicon	Yamato Scientific
Refrigerator switch	2-01-027-0005	HLN215A	Fujisoku
Trouble lamp	2-09-006-0001	BN-5701	Satoh Parts
Earth leakage breaker	2-06-000-0019	FG32R/15-30MA15A	Fuji Denki
Drain hose	3-04-001-6022	ϕ 6 \times ϕ 18	Yamato Scientific
Drain tap	CA300-40200	Polyacetal	Yamato Scientific
Compressor	3-01-006-0005	C-2SN350LOR HFC23/R404A	Yamato Scientific
Oil separator	3-19-000-0003	OUB-1-S-3	Yamato Scientific
Pressure switch	3-18-000-0006	VFP-F	Fuji Kohki
Condenser	CA800-30140		Yamato Scientific
Fan motor	3-01-006-0006	SE4-CO41NP	Sanyo
Charge bulb	3-25-001-0002		Yamato Scientific
Dryer	3-20-003-6002	KC-10432	Meikoh Kiki
Terminal block	LT00035673	MKH-250ABC-5P	Terminal
Relay	2-05-008-0002	AP3124K	Matsushita
Relay	2-05-000-0035	LY2F AC100V	Omron
Time relay	2-05-000-0053	ADX11134	Matsushita
Socket	2-05-008-1004	AP3822K	Matsushita
Vacuum hose	3-04-001-6022	ϕ 6 \times ϕ 18	Yamato Scientific

List of Dangerous Substances



Never use explosive substances, flammable substances and substances that include explosive or flammable ingredients in this unit.

EXPLOSIVE

EXPLOSIVE:	Ethylene glycol dinitrate (nitro glycol), Glycerin trinitrate (nitroglycerine), Cellulose nitrate (nitrocellulose), and other explosive nitrate esters
	Trinitrobenzene, Trinitrotoluene, Trinitrophenol (picric acid), and other explosive nitro compounds
	Acetyl hidroperoxide (peracetic acid), Methyl ethyl ketone peroxide, Benzyl peroxide, and other organic peroxides

FLAMMABLE

IGNITING:	Lithium (metal), Potassium (metal), Sodium (metal), Yellow phosphorus, Phosphorus sulfide, Red phosphorus, Celluloid compounds, Calcium carbide, Lime phosphate, Magnesium (powder), Aluminum (powder), Powder of metals other than magnesium and aluminum, Sodium hydrosulfite
OXIDIZING:	Potassium chlorate, Sodium chlorate, Ammonium chlorate, and other chlorate
	Potassium perchlorate, Sodium perchlorate, Ammonium perchlorate, and other perchlorate
	Potassium peroxide, Sodium peroxide, Barium peroxide, and other inorganic peroxide
	Potassium nitrate, Sodium nitrate, Ammonium nitrate, and other nitrate
	Sodium chlorite and other chlorites
Calcium hypochlorite and other hypochlorites	
INFLAMMABLE LIQUID:	Ethyl ether, Gasoline, Acetaldehyde, Propylene chloride, Carbon disulfide, and other flammable substances having a flash point of lower than -30°C
	Normal hexane, ethylene oxide, acetone, benzene, methyl ethyl ketone, and other flammable substances having a flash point of -30°C or higher but lower than 0°C
	Methanol, Ethanol, Xylene, Pentyl acetate (amyl acetate), and other flammable substances having a flash point of 0°C or higher but lower than 30°C
	Kerosene, Light oil (gas oil), Oil of turpentine, Isopentyl alcohol (isoamyl alcohol), Acetic acid, and other flammable substances having a flash point of 30°C or higher but lower than 65°C
FLAMMABLE GAS:	Hydrogen, Acetylene, Ethylene, Methane, Propane, Butane, and other flammable substances which assume a gaseous state at 15°C and 1 atm

(Source: Appendix Table 1 of Article 6 of the Industrial Safety and Health Order in Japan)

Responsibility

Please follow the instructions in this document when using this unit. Yamato Scientific has no responsibility for the accidents or breakdown of device if it is used with a failure to comply. Never conduct what this document forbids. Unexpected accidents or breakdown may result in.

Note

- ◆ The contents of this document may be changed in future without notice.
- ◆ Any books with missing pages or disorderly binding may be replaced.

Instruction Manual for

Cold Trap

Model CA300/800

First Edition Feb. 21, 2002

Revised May. 11, 2012

Yamato Scientific Co., Ltd.

2-1-6 Nihonbashi Honcho, Chuo-ku,
Tokyo, 103-8432, Japan

<http://www.yamato-net.co.jp>